


EXTREMIFIX™

Cannulated Screw System | Mini & Small

Product Information

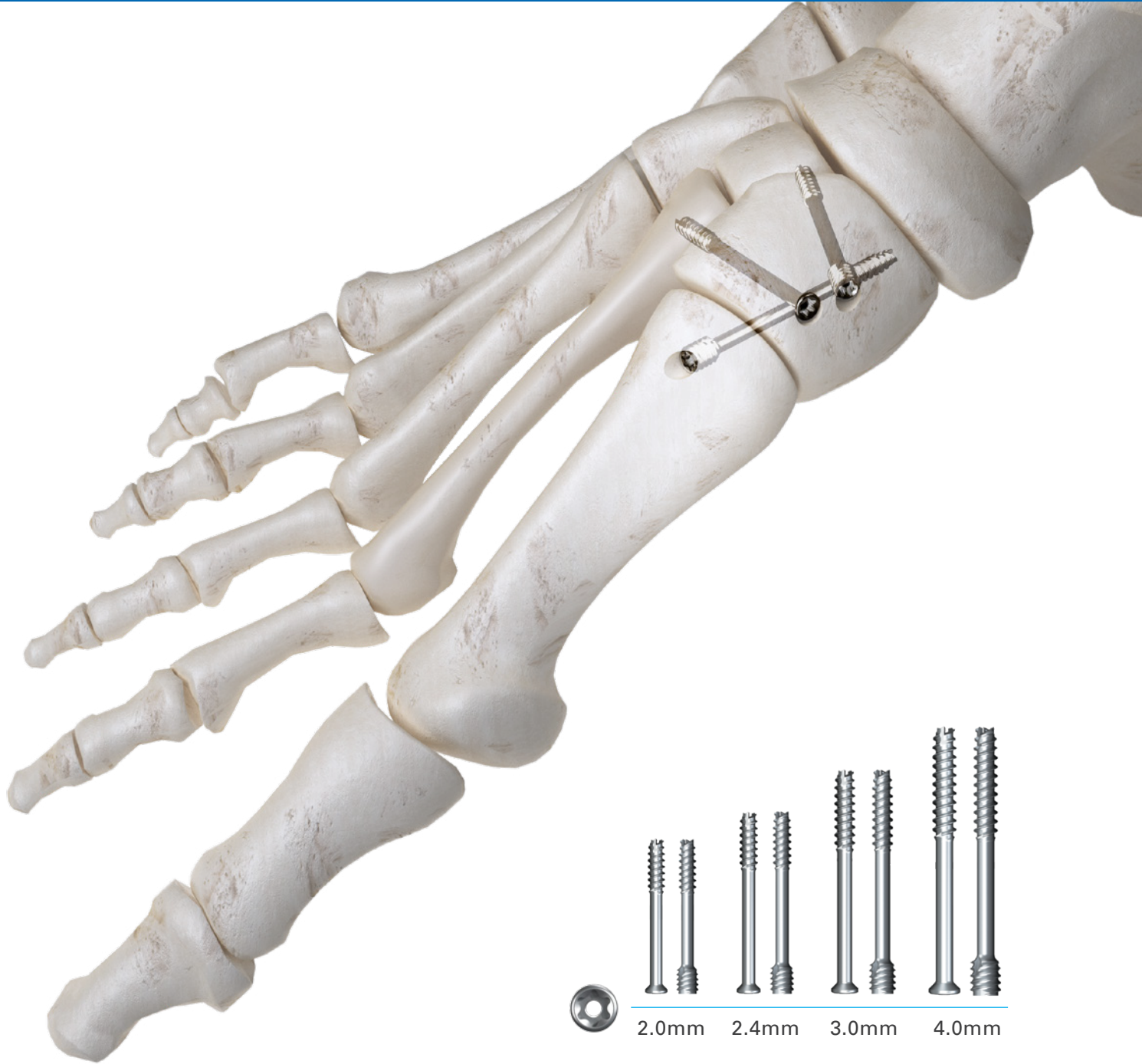

Rethinking Possibilities, Reshaping Lives

ExtremiFix Mini & Small

The ExtremiFix Mini & Small Cannulated Screw System offers multiple screw and instrument options and sizes for fracture fixation, fusion, and osteotomy procedures of the upper and lower extremities. The system offering includes Ø2.0mm, Ø2.4mm, Ø3.0mm, and Ø4.0mm partially-threaded, headed and headless screws in various lengths.

Headed Screw Design

Headless Screw Design

Low-Profile Head
Minimizes soft tissue irritation

Conical Head
Minimizes hardware prominence after countersinking and creates concentric compression forces



Hexalobe Drive Recess T7 and T10
Higher torque resistance

Reverse Cutting Flutes
Easy removal

Double Lead
(Ø3.0mm & Ø4.0mm only)
Faster advancement

Sharp Threads
Bony engagement,
purchase, stability, fixation

Self-Drilling and Self-Tapping
Easy insertion

Headless Screw Compression
Differential pitch between the cancellous, greater-pitched, leading threads and the cortical, finer-pitched, trailing threads

Minimize Protrusion
Allows for flush placement within the bone



Washers

The ExtremiFix Mini & Small Cannulated Screw System includes washers for use with the Ø2.0mm, Ø2.4mm, Ø3.0mm, and Ø4.0mm headed screws. Washers may be used to provide additional contact area in patients with soft bone and are laser-marked to indicate which side must be in contact with the screw head.



Refer to the provided instructions for use for the complete indications, contraindications, warnings, and instructions for use.



3885 Arapaho Rd.
Addison, TX 75001
Customer Service: 800.456.7779
Outside the U.S.: 001.972.677.4600
Fax: 800.390.2620
Fax Outside the U.S.: 001.972.677.4709
E-mail: customer.service@osteomed.com
www.osteomed.com